


iM3 - The Veterinary Dental Company



iM3 elite LED Dental Unit

The elite dental system features advanced LED light in the high speed handpiece, automatic handpiece activators and the only fully autoclavable air driven dental suction in the industry.

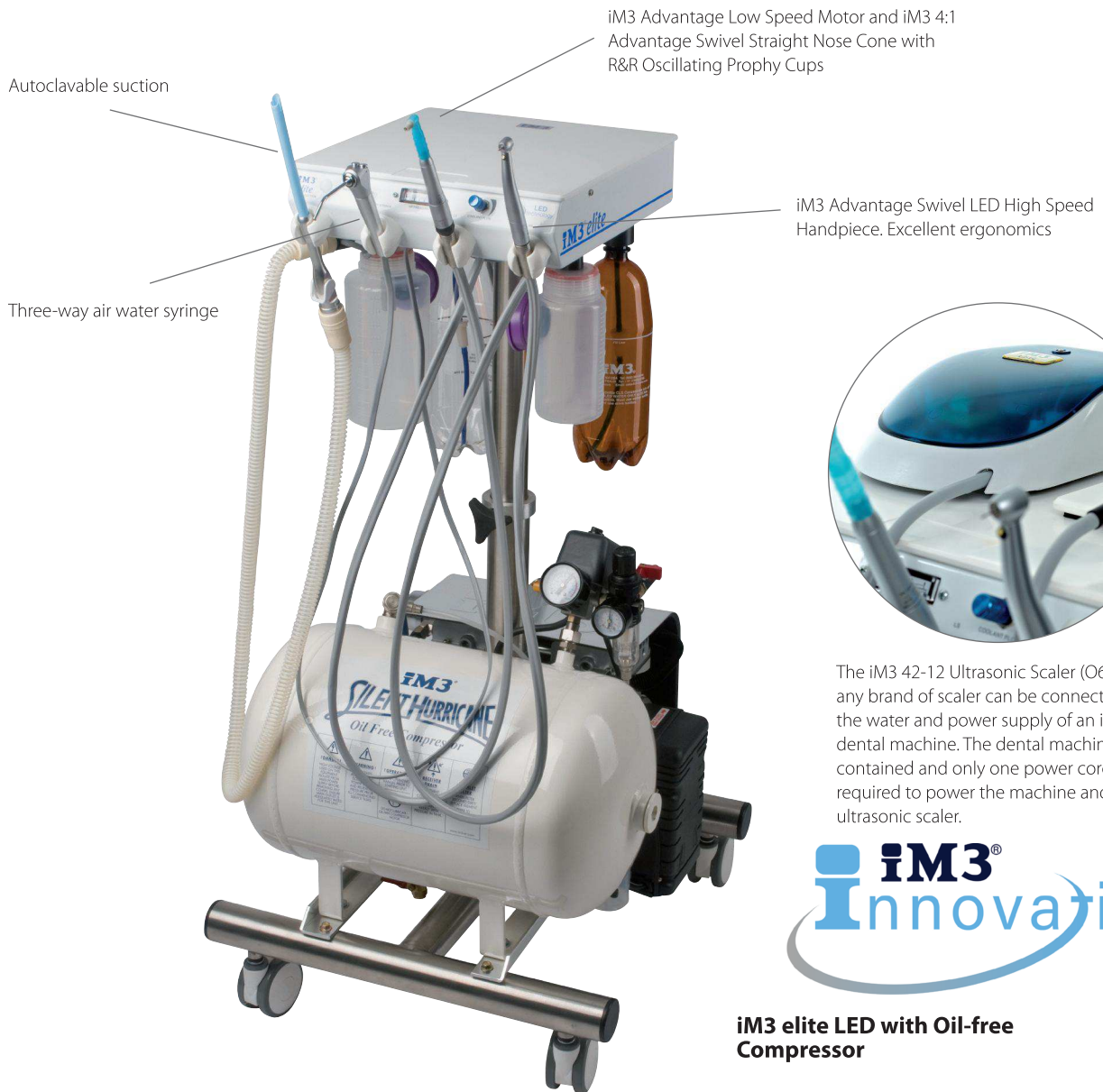

Autoclavable suction

Three-way air water syringe

iM3 Advantage Low Speed Motor and iM3 4:1 Advantage Swivel Straight Nose Cone with R&R Oscillating Prophy Cups

iM3 Advantage Swivel LED High Speed Handpiece. Excellent ergonomics



The iM3 42-12 Ultrasonic Scaler (O6979) or any brand of scaler can be connected to the water and power supply of an iM3 dental machine. The dental machine is self contained and only one power cord is required to power the machine and ultrasonic scaler.



iM3 elite LED with Oil-free Compressor

Cat. No
U3500

Don't forget to add iM3 Straw

Cat. No
PS001

iM3 - The Veterinary Dental Company



Features and Benefits of the iM3 elite Dental Unit

iM3 has developed an air-driven dental suction unit that is fully autoclavable. The suction tips, handpiece and smooth bore suction tubing can be easily removed and autoclaved.

The suction waste is collected in an easily accessible one litre container with a replaceable filter. A separate foot pedal operates the suction, allowing its use in conjunction with the high speed drill or three-way air, water, and mist syringe. Because it's fully autoclavable, it can also be used in surgery.

The quality of the handpieces and iM3's unique features ensure the iM3 elite Dental Unit is regarded as the number one in the industry.

The iM3 elite comes with:

- Suction
- Three-way syringe (air, water and mist)
- Advantage Low Speed handpiece
- Advantage Swivel LED High Speed Handpiece

The iM3 designed magnetic bur holder prevents burs from falling on the ground. The burs are held in place right where you need them. The Corian® surface is impervious to bacteria and easy to clean unlike fibreglass and plastics.

The dental machine chassis is made from aluminium and hospital-grade powder coating. The unit's control box and components are covered by a full five-year warranty. The iM3 elite, with its built-in IEC power board and coolant outlet, allows for the easy integration of an ultrasonic scaler.

The iM3 elite dental unit can be supplied with a height adjustable (50 cm) long reach (150 cm) pneumatic wall arm. The arm allows the machine to be stored against the wall and extended to the end of the wet table for use.

The air line and electricity for the dental machine are contained in the arm, removing unnecessary cables from the ground.

All iM3 dental units are fitted with an iM3 Advantage Low Speed Motor and iM3 4:1 Advantage Swivel Straight Nose Cone for polishing. An array of low speed attachments for rabbit dentistry, orthopaedics and restorative dentistry are available.

Not all automatic handpiece switches are the same! iM3 has developed a pneumatic switch with no electric parts. The switch removes the need to manually select the high or low speed handpiece. If it ever needs to be replaced, it can be done in seconds by the user, which gives peace of mind.



iM3 elite LED with Oil-free Compressor

Cat. No
U3500

iM3 elite LED SW

(Slow speed with water cooling, not available in all countries)

U3500SW

Don't forget to add iM3 Straw

Cat. No
PS001