

iM3 - The Veterinary Dental Company

Maintenance instructions conveniently printed on the back of the dental machine

High speed handpiece for fast sectioning of teeth

Low speed polisher

Three-way air water syringe for rinsing the oral cavity

Adjustable stainless steel stand



iM3 oil-free compressor with large tank so it doesn't overheat

iM3 GS Deluxe START Dental Unit

The GS START is a professional dental system with endless options for the future. The START was designed as a way of getting "STARTed" in veterinary dentistry with the opportunity to upgrade the machine over time as your dental knowledge and dental treatments increase. Rather than buy a cheap dental unit or electric micro motor that won't last, buy a START and let it grow with your practice!

The START can evolve from having a simple one bottle water system to a more advanced two-bottle system with CLS at your fingertips. Then upgrade to an Advantage Swivel LED High Speed Handpiece with swivel and triple coolant spray.

The GS Deluxe START

- Air driven high speed and low speed handpieces
- 400,000 rpm Advantage highspeed handpiece with coolant control
- Low speed polisher
- Three-way syringe for air, water and mist
- Oil-free Silent Hurricane low-maintenance compressor
- Upgrade to GS Deluxe or Deluxe LED at any time

The START is not a cheap dental unit. It is a high quality iM3 built dental system with the same quality valves and components used in all iM3 machines. The START is designed to offer all the features needed on an entry level dental machine at the same time, providing the quality and longevity you have come to expect from iM3.




iM3 GS Deluxe START with Oil-free Compressor

Cat. No
U3010

Don't forget to add iM3 Straw

Cat. No
PS001

iM3 - The Veterinary Dental Company

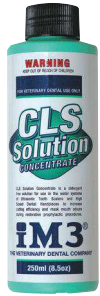

The GS Deluxe dental machine has a two bottle coolant system. The clear bottle is filled with distilled water, while the brown bottle protects the UV light-sensitive CLS solution. CLS will help minimise cross contamination when used between patients. The spearmint smell will mask mouth odours.

Three-way syringe - ideal for washing away prophy paste and debris

iM3 GS Deluxe Dental Unit

The GS Deluxe includes a high-quality Advantage High Speed handpiece (400,000 rpm) with triple-coolant spray for fast sectioning of teeth; a three-way syringe with the ability to select water, CLS, air, or a mist; and a low speed handpiece for polishing using iM3 Rock & Roll oscillating prophy angles that won't wind up muzzle hair.

The three-way (air, water and mist) syringe is a useful instrument for washing away debris from the mouth. The tip is removable and autoclavable and can be rotated 360 degrees. The use of CLS in the dental machine allows for easy application of chlorhexidine into the mouth during dental procedures. Stainless steel instrument tray with silicon mat provides added storage and the ideal location for your iM3 dental instruments.

The iM3 42-12 Ultrasonic Scaler or any brand of scaler can be connected to the water and power supply of the GS Deluxe dental machine. The dental machine is self-contained and only one power cord is required to power the machine and ultrasonic scaler.



Low speed polisher with oscillating prophy angles that won't catch muzzle hair

iM3 Advantage High Speed handpiece for fast sectioning of teeth and bone removal

iM3 five-year warranty with unlimited factory direct technical support



iM3 GS Deluxe (non LED) with Oil-free Compressor

Cat. No
U3012

Don't forget to add iM3 Straw

Cat. No
PS001